




3-Furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

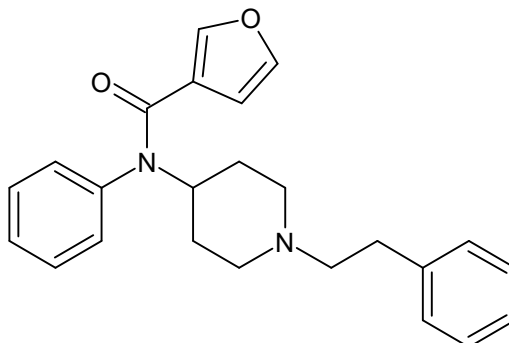

1. GENERAL INFORMATION

IUPAC Name: N-phenyl-N-[1-(2-phenylethyl)piperidin-4-yl]furan-3-carboxamide

CAS#: NA

Synonyms: NA

Source: DEA Reference Material Collection

Appearance: White powder

UV_{max}(nm): Not determined

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C ₂₄ H ₂₆ N ₂ O ₂	374.48	Not Determined
HCl	C ₂₄ H ₂₆ N ₂ O ₂ HCl	410.94	Not Determined



3-Furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~15 mg/mL in Methanol - d_4 containing TMS for 0 ppm reference and dimethylfumarate as quantitative internal standard.

Instrument: 400 MHz NMR spectrometer
Parameters: Spectral width: at least containing -3 ppm through 13 ppm
Pulse angle: 90°
Delay between pulses: 45 seconds

^1H NMR: 3-Furanyl fentanyl HCl; Lot # ALB-288-6; Methanol - d_4 ; 400MHz

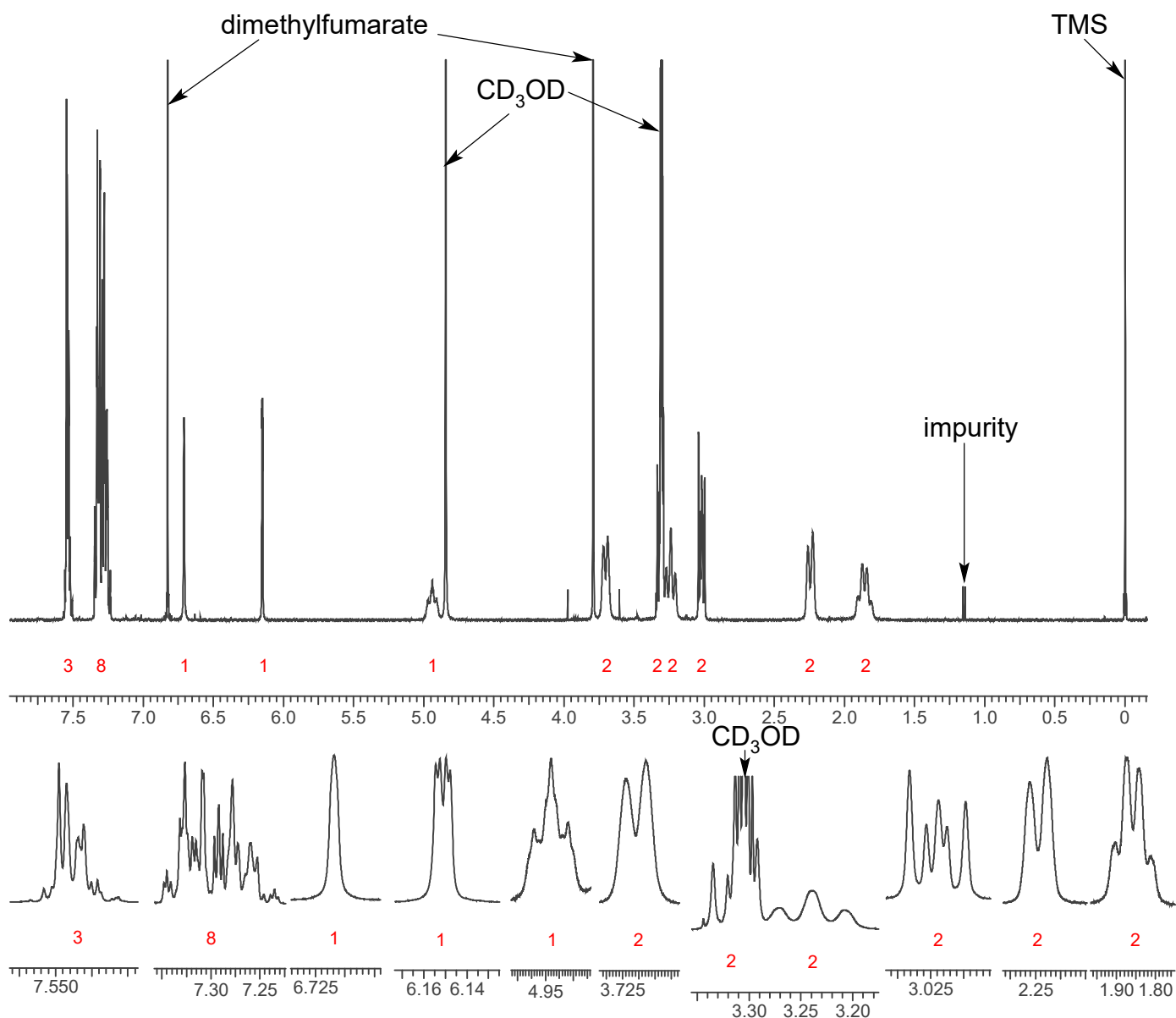



3-Furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

Sample Preparation: Dilute analyte ~4 mg/mL in CHCl₃

Instrument: Agilent gas chromatograph operated in split mode with MS detector
Column: DB-5 MS (or equivalent); 15m x 0.25 mm x 0.25 μm
Carrier Gas: Helium at 1.5 mL/min
Temperatures: Injector: 280°C MSD transfer line: 280°C
MS Source: 250°C MS Quad: 150°C
Oven program:
 1) 100°C initial temperature for 1.0 min
 2) Ramp to 280°C at 12 °C/min
 3) Hold final temperature for 9.0 min
Injection Parameters: Split Ratio = 25:1, 1 μL injected
MS Parameters: Mass scan range: 30-550 amu Threshold: 150
Tune file: stune.u Acquisition mode: scan
Retention Time: 16.65 min

EI Mass Spectrum: 3-Furanyl fentanyl HCl; Lot # ALB-288-6

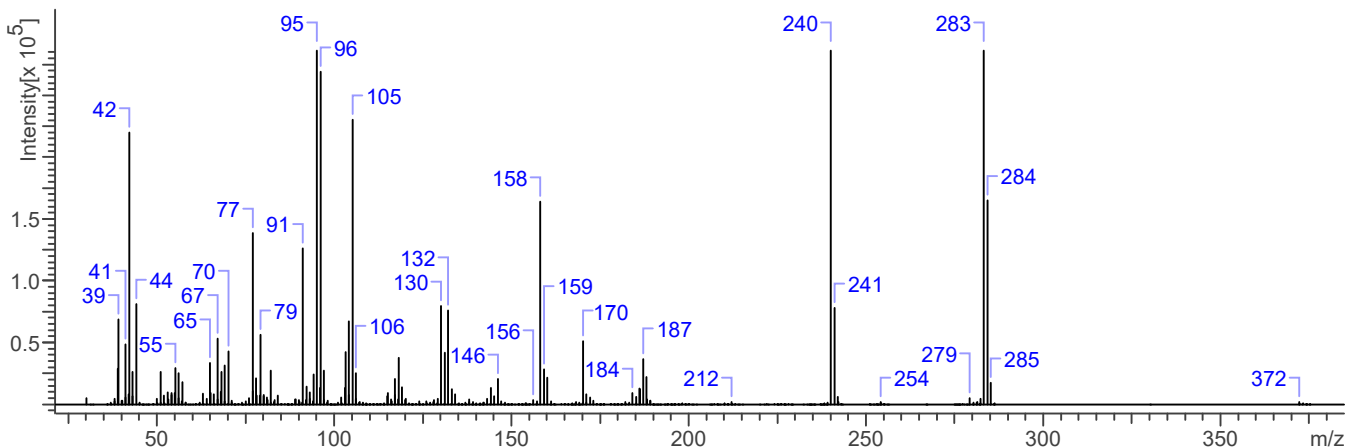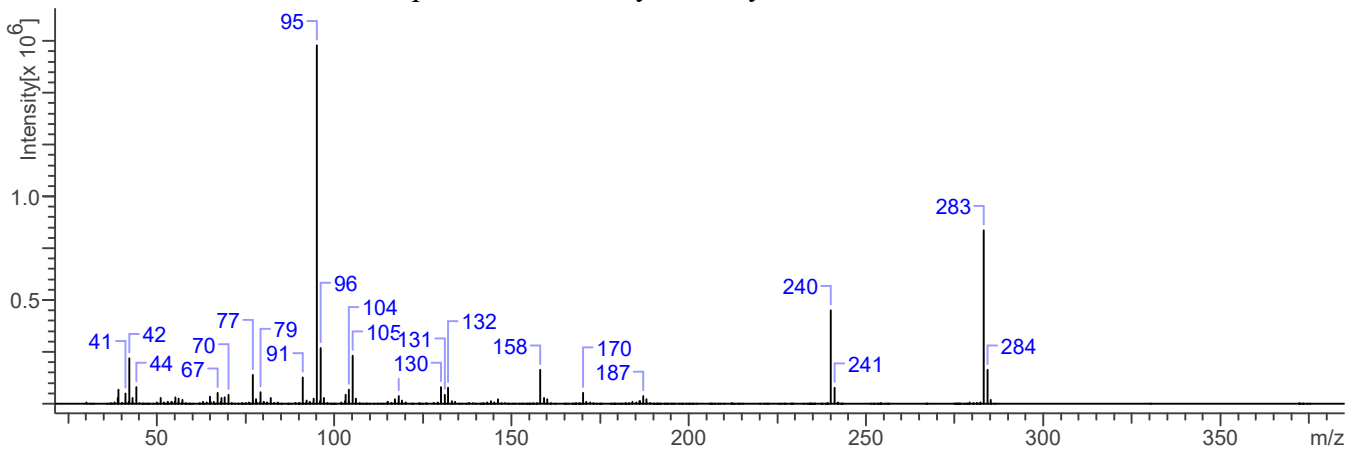



3-Furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with diamond ATR attachment (1 bounce)

Scan Parameters:
Number of scans: 32
Number of background scans: 32
Resolution: 4 cm⁻¹
Sample gain: 4
Aperture: 80

FTIR ATR (Diamond 1 Bounce): 3-Furanyl fentanyl HCl; Lot # ALB-288-6

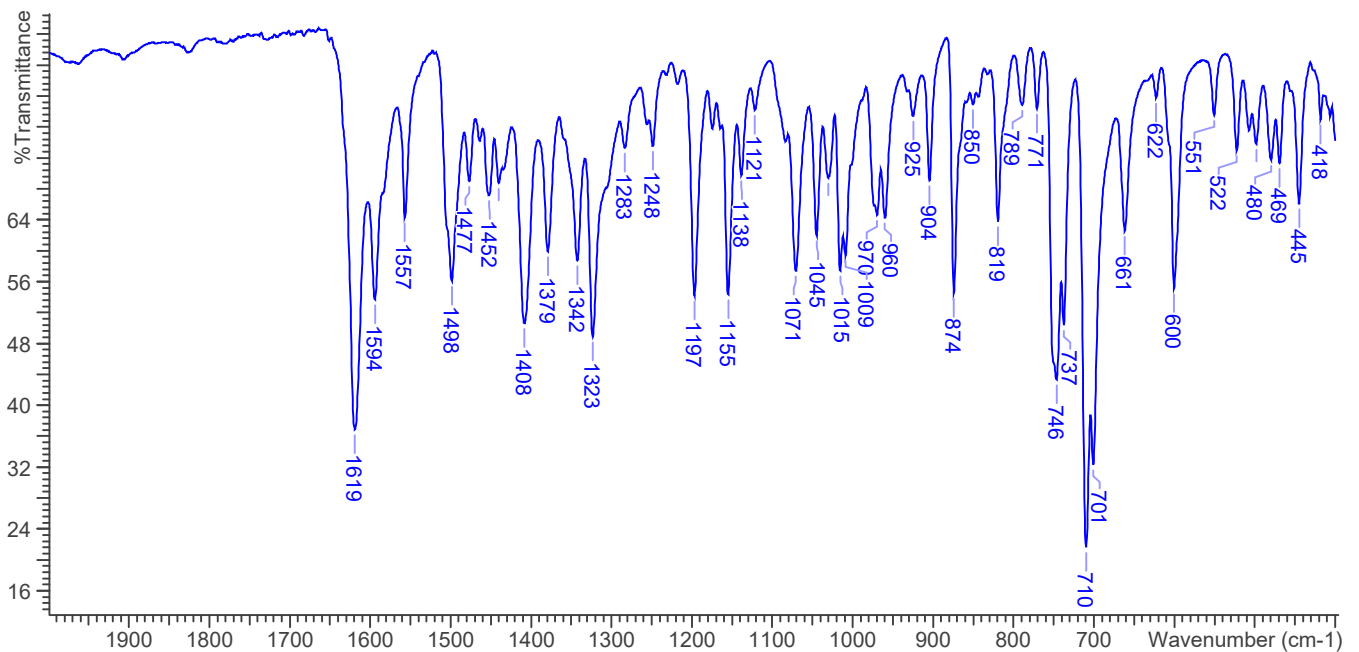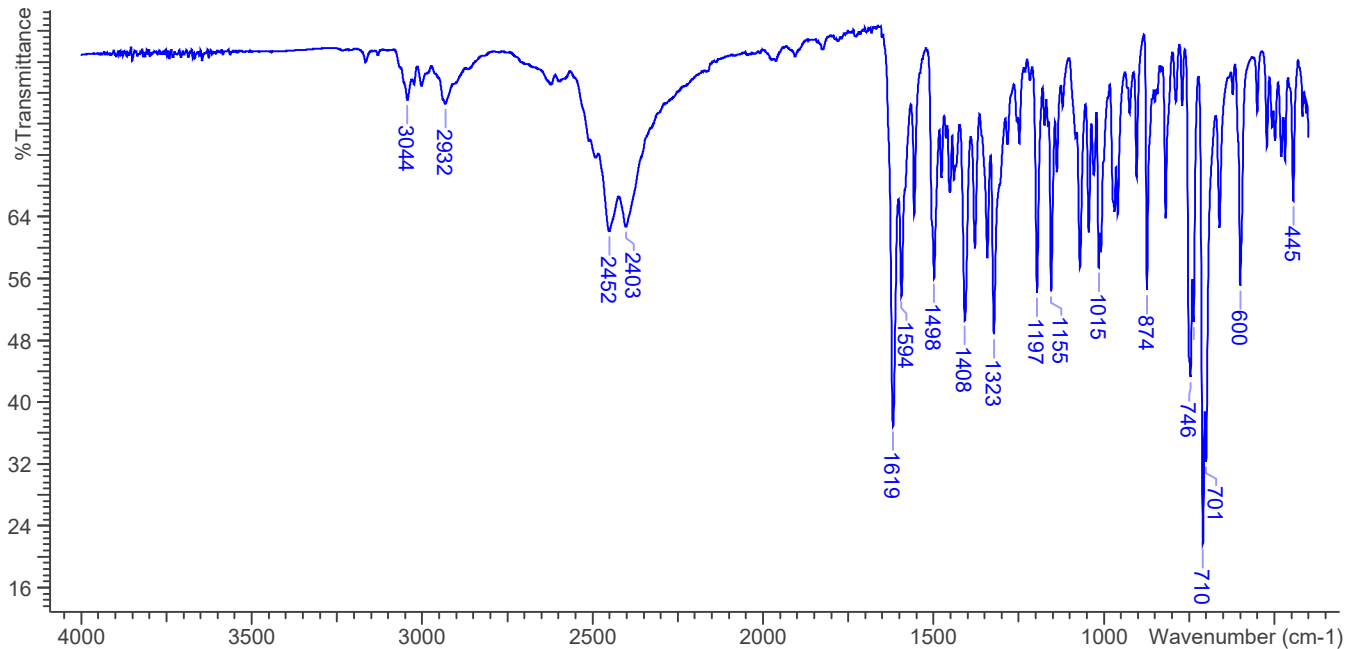



3-Furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



4. ADDITIONAL RESOURCES

No additional resources as of 09/2016