




Cyclopentanoyl Fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

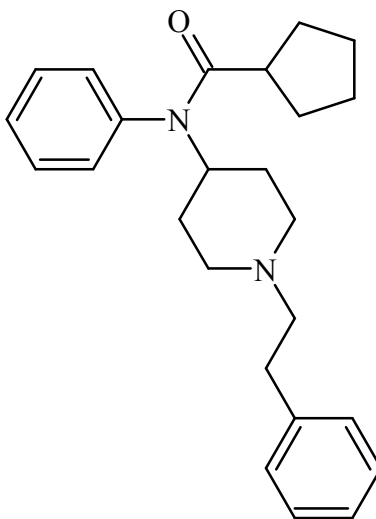

1. GENERAL INFORMATION

IUPAC Name:	N-phenyl-N-[1-(2-phenylethyl)piperidin-4-yl]cyclopentanecarboxamide
CAS#:	N/A
Synonyms:	N/A
Source:	DEA Reference Material Collection
Appearance:	Off-White powder
UV_{max}(nm):	Not Determined

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C ₂₅ H ₃₂ N ₂ O	376.53	Not Determined
HCl	C ₂₅ H ₃₂ N ₂ O HCl	413.00	Not Determined



Cyclopentanoyl Fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~10 mg/mL in CD₃OD containing TMS for 0 ppm reference and dimethylfumarate as quantitative internal standard.


Instrument: 400 MHz NMR spectrometer

Parameters: Spectral width: at least containing -3.0 ppm through 13.0 ppm

Pulse angle: 90°

Delay between pulses: 45 seconds

¹HNMR: Cyclopentanoyl fentanyl HCl; Lot ALB-297-8-2; CD₃OD; 400MHz

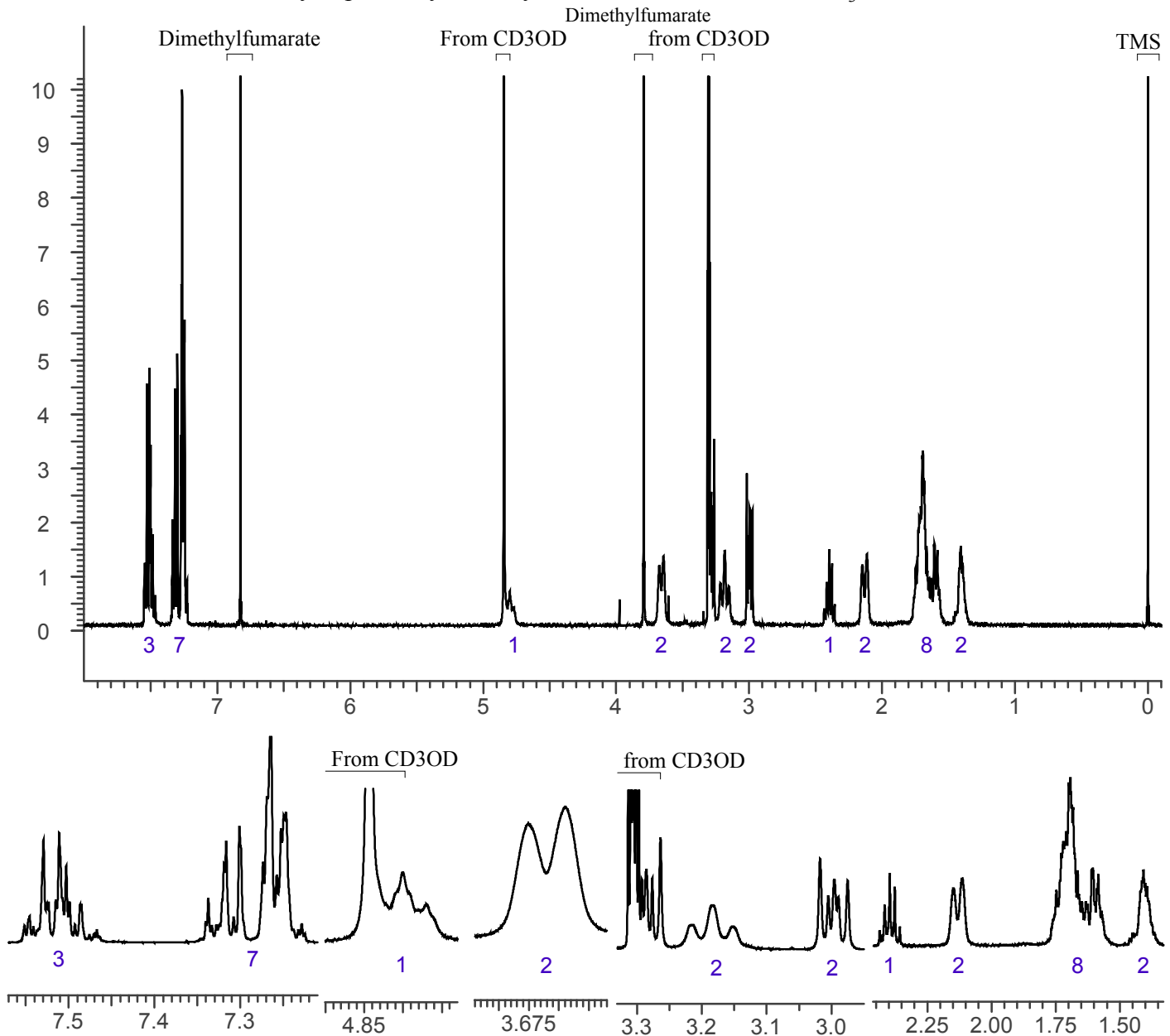



Cyclopentanoyl Fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

Sample Preparation: Dilute analyte ~4 mg/mL in chloroform.

Instrument: Agilent gas chromatograph operated in split mode with MS detector

Column: HP-5 MS (or equivalent); 30m x 0.25 mm x 0.25 μ m

Carrier Gas: Helium at 1.5 mL/min

Temperatures: Injector: 280°C

MSD transfer line: 280°C

MS Source: 230°C

MS Quad: 150°C

Oven program:

1) 100°C initial temperature for 1.0 min

2) Ramp to 280°C at 12 °C/min

3) Hold final temperature for 9.0 min

Injection Parameters: Split Ratio = 25:1, 1 μ L injected

MS Parameters: Mass scan range: 30-550 amu



Threshold: 100

Tune file: stune.u

Acquisition mode: scan

Retention Time: 20.3 min

EI Mass Spectrum: Cyclopentanoyl Fentanyl HCl; Lot ALB-297-8-2

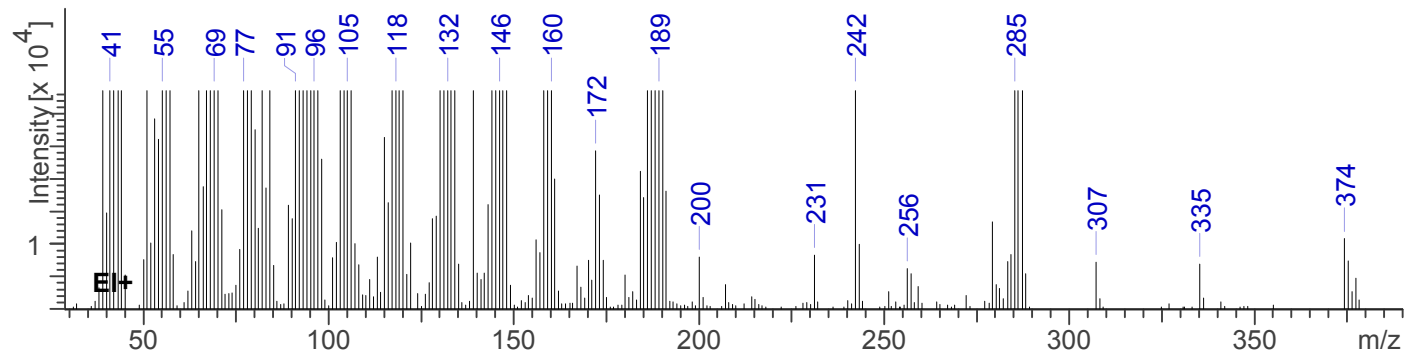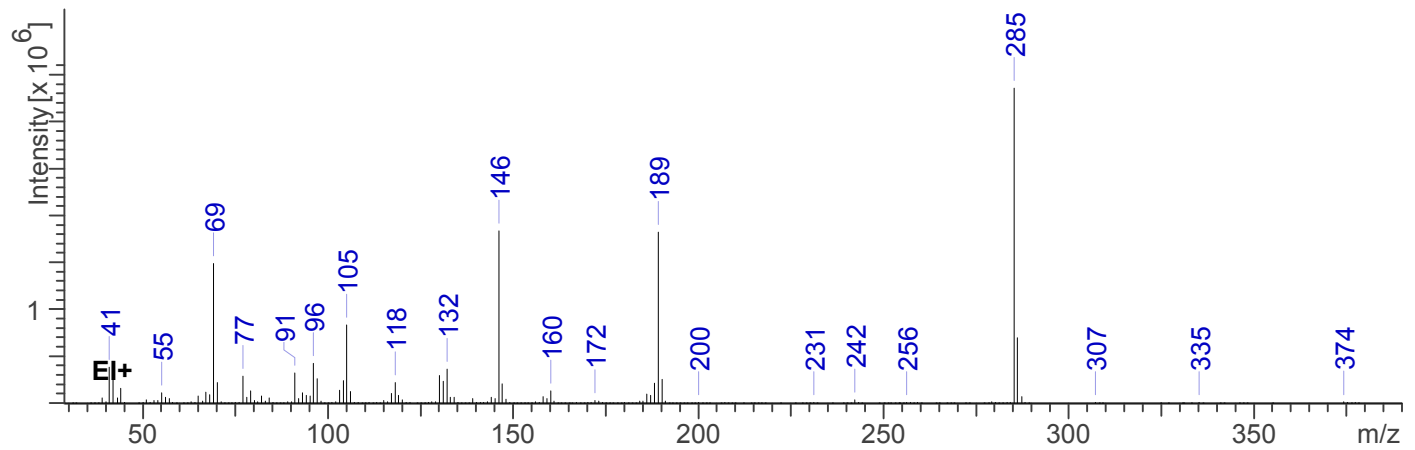



Cyclopentanoyl Fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with Golden Gate diamond ATR attachment (3 bounce)

Scan Parameters:
Number of scans: 16
Number of background scans: 16
Resolution: 4 cm^{-1}
Sample gain: 8
Aperture: 150

FTIR (Golden Gate with ATR, 3 Bounce): Cyclopentanoyl Fentanyl HCl; Lot ALB-297-8-2

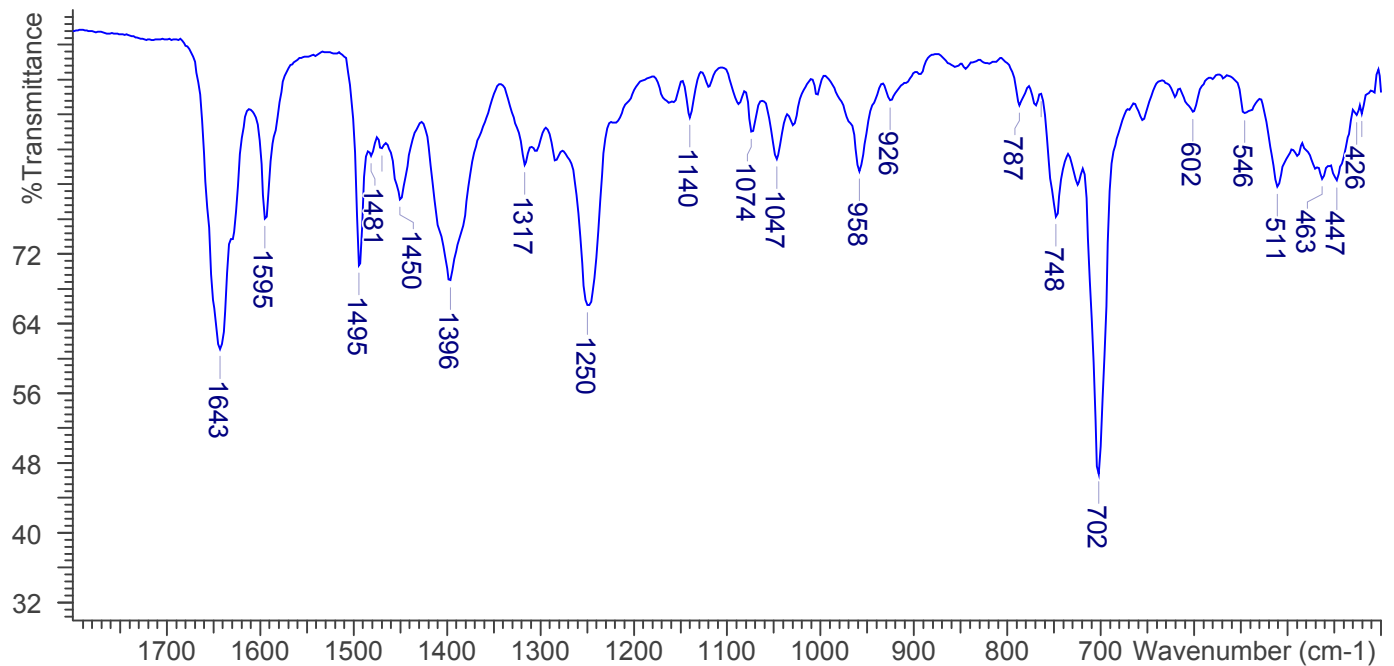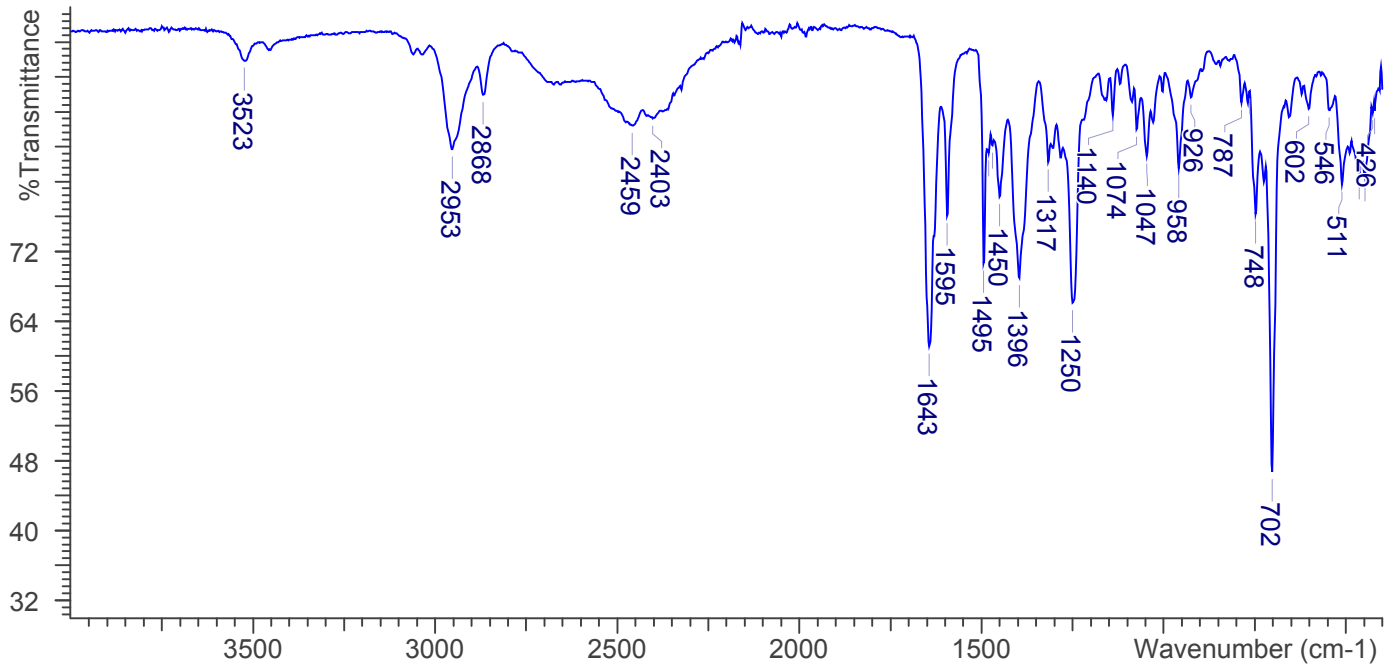



Cyclopentanoyl Fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



4. ADDITIONAL RESOURCES

No Literature as of 08/2017