




alpha-benzyl-N-methylphenethylamine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

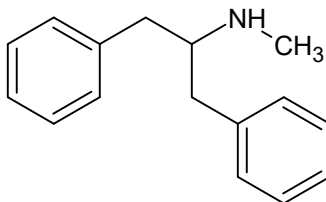

1. GENERAL INFORMATION

<i>IUPAC Name:</i>	<i>N</i> -methyl-1,3-diphenylpropan-2-amine
<i>CAS#:</i>	N/A
<i>Synonyms:</i>	BNMPA, B-compound
<i>Source:</i>	DEA Reference Material Collection
<i>Appearance:</i>	White powder
<i>UV_{max}(nm):</i>	Not determined

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C ₁₆ H ₁₉ N	225.33	Not Determined
HCl	C ₁₆ H ₁₉ N HCl	261.79	194.83



alpha-benzyl-*N*-methylphenethylamine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~15 mg/mL in DMSO- d_6 containing TMS for 0 ppm reference and maleic acid as quantitative internal standard.

Instrument: 400 MHz NMR spectrometer
Parameters: Spectral width: at least containing -3 ppm through 13 ppm
Pulse angle: 90°
Delay between pulses: 45 seconds

^1H NMR: *alpha*-benzyl-*N*-methylphenethylamine HCl; Lot# 2TDM-208-01; DMSO- d_6 ; 400MHz

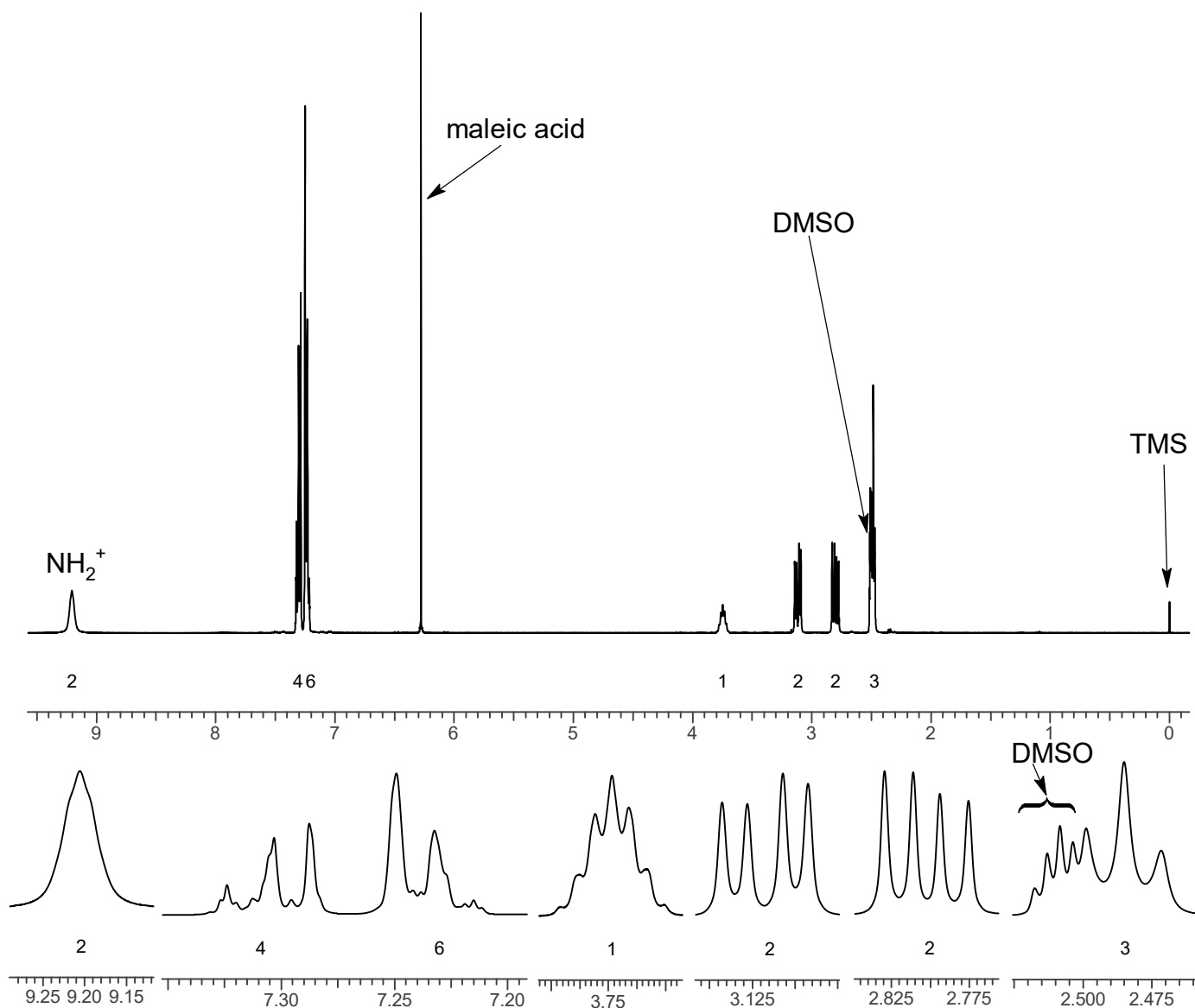



alpha-benzyl-*N*-methylphenethylamine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

Sample Preparation: Dilute analyte ~4 mg/mL in CHCl₃, base extracted

Instrument: Agilent gas chromatograph operated in split mode with MS detector
Column: HP-5 MS (or equivalent); 30m x 0.25 mm x 0.25 μm
Carrier Gas: Helium at 1.5 mL/min
Temperatures: Injector: 280°C MSD transfer line: 280°C
MS Source: 230°C MS Quad: 150°C
Oven program:
1) 100°C initial temperature for 1.0 min
2) Ramp to 280°C at 12 °C/min
3) Hold final temperature for 9.0 min
Injection Parameters: Split Ratio = 25:1, 1 μL injected
MS Parameters: Mass scan range: 30-550 amu Threshold: 250
Tune file: stune.u Acquisition mode: scan
Retention Time: 10.05 min

EI Mass Spectrum: *alpha*-benzyl-*N*-methylphenethylamine HCl; Lot# 2TDM-208-01, base extracted

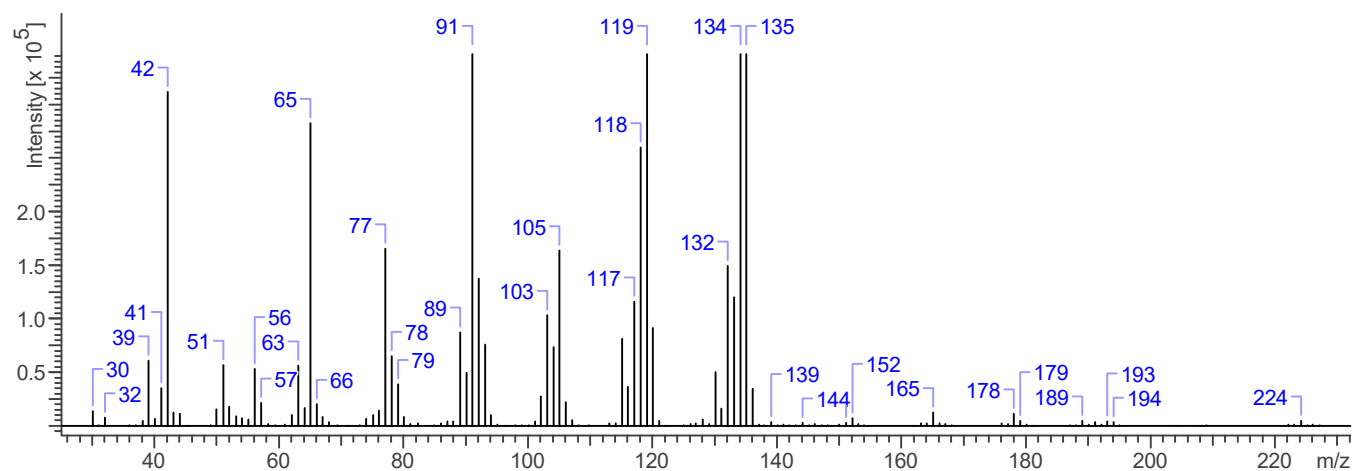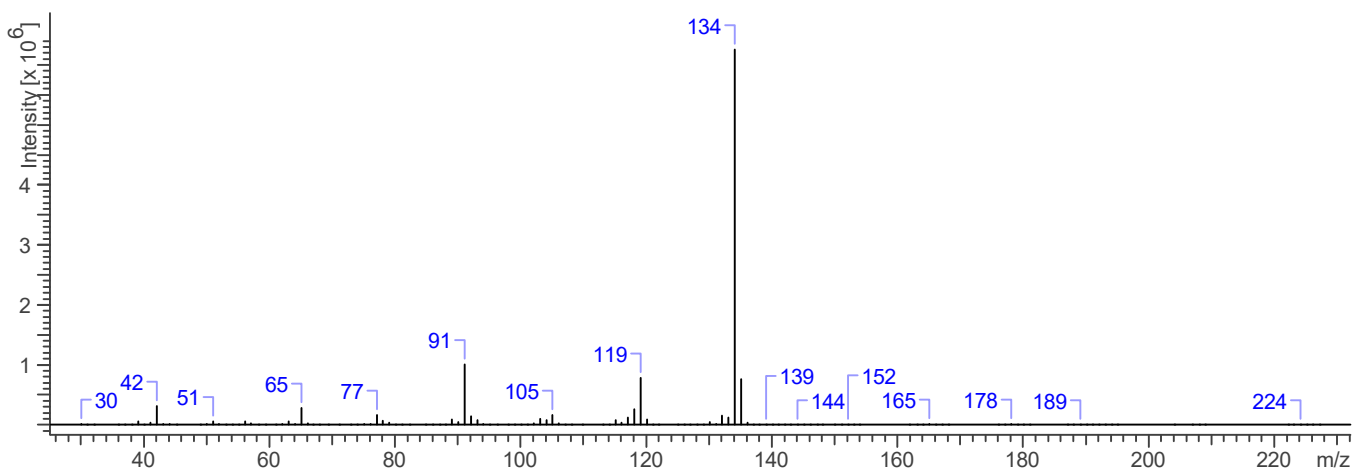



alpha-benzyl-*N*-methylphenethylamine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with diamond ATR attachment (1 bounce)

Scan Parameters:
Number of scans: 32
Number of background scans: 32
Resolution: 4 cm⁻¹
Sample gain: 1
Aperture: 150

FTIR ATR (Diamond 1 Bounce): *alpha*-benzyl-*N*-methylphenethylamine HCl; Lot# 2TDM-208-01

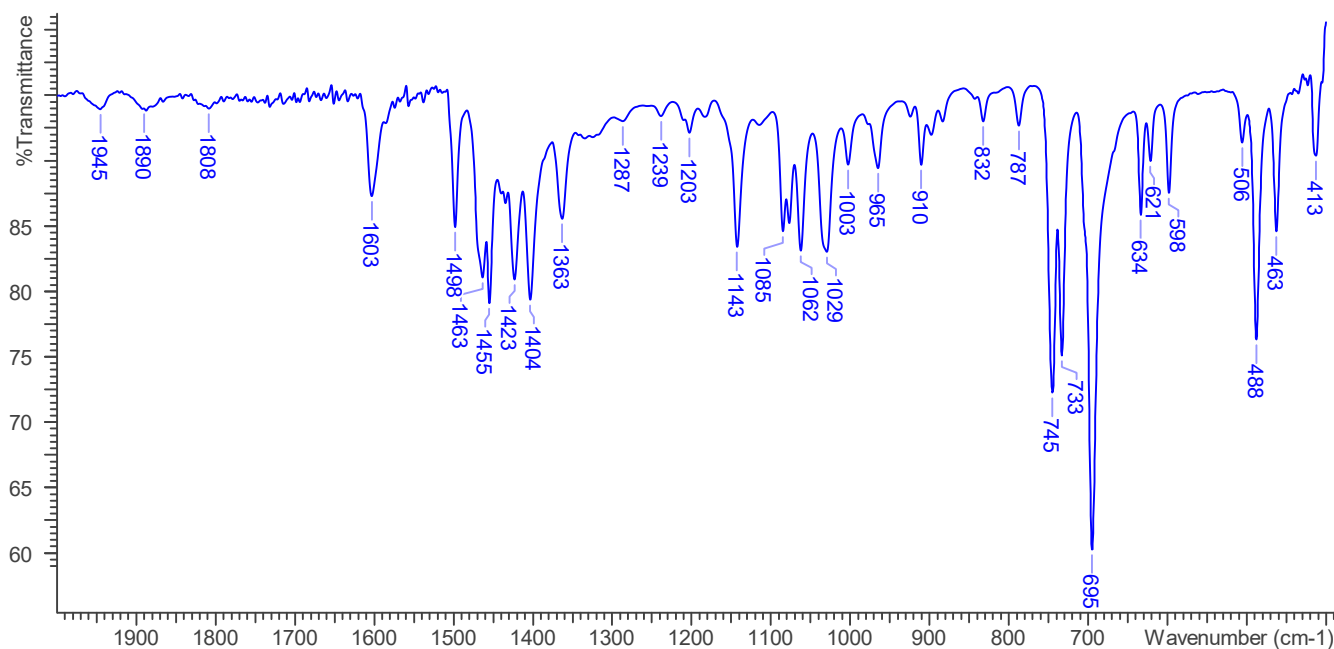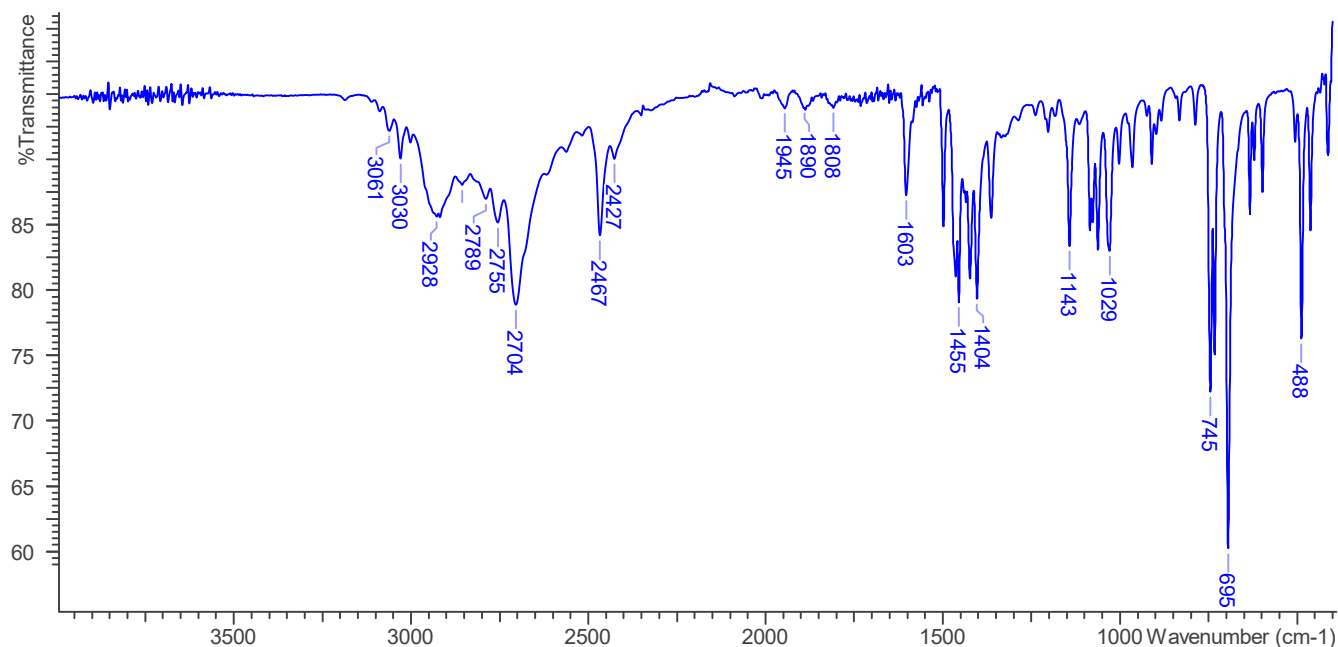



alpha-benzyl-N-methylphenethylamine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



4. ADDITIONAL RESOURCES

1. Barron, Robert P., et al, Identification of Impurities in Illicit Methamphetamine Samples. J Assoc Off Anal Chem. 1974, 57(5):1147-58.
2. Kram, T. C.; Kruegel, A. V., The Identification of Impurities in Illicit Methamphetamine Exhibits by Gas Chromatography/Mass Spectrometry and Nuclear Magnetic Resonance Spectroscopy. J Forensic Sciences 1977, 22(1):40-52.