




para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

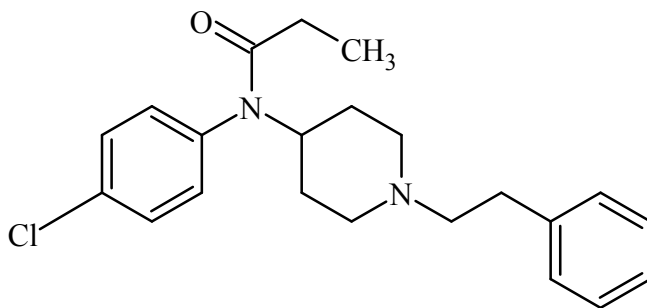

1. GENERAL INFORMATION

IUPAC Name:	N-(4-chlorophenyl)-N-(1-phenethylpiperidin-4-yl)propanamide
CAS#:	NA
Synonyms:	N-(4-chlorophenyl)-N-[1-(2-phenylethyl)piperidin-4-yl]propanamide, N-(4-chlorophenyl)-N-[1-(2-phenylethyl)-4-piperidinyl]-propanamide, 4-Chlorofentanyl, p-Chlorofentanyl
Source:	DEA Reference Material Collection
Appearance:	white powder
UV_{max}(nm):	NA

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C ₂₂ H ₂₇ ClN ₂ O	370.92	NA
HCl	C ₂₂ H ₂₇ ClN ₂ O · HCl	407.38	NA



para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~11 mg/mL in CD₃OD containing TMS for 0 ppm reference and 1,4-BTMSB-d₄ as quantitative internal standard.

Instrument: 400 MHz NMR spectrometer
Parameters: Spectral width: at least containing -3 ppm through 13 ppm
Pulse angle: 90°
Delay between pulses: 45 seconds

¹HNMR: para-Chlorofentanyl HCl; Lot# 0560132-2

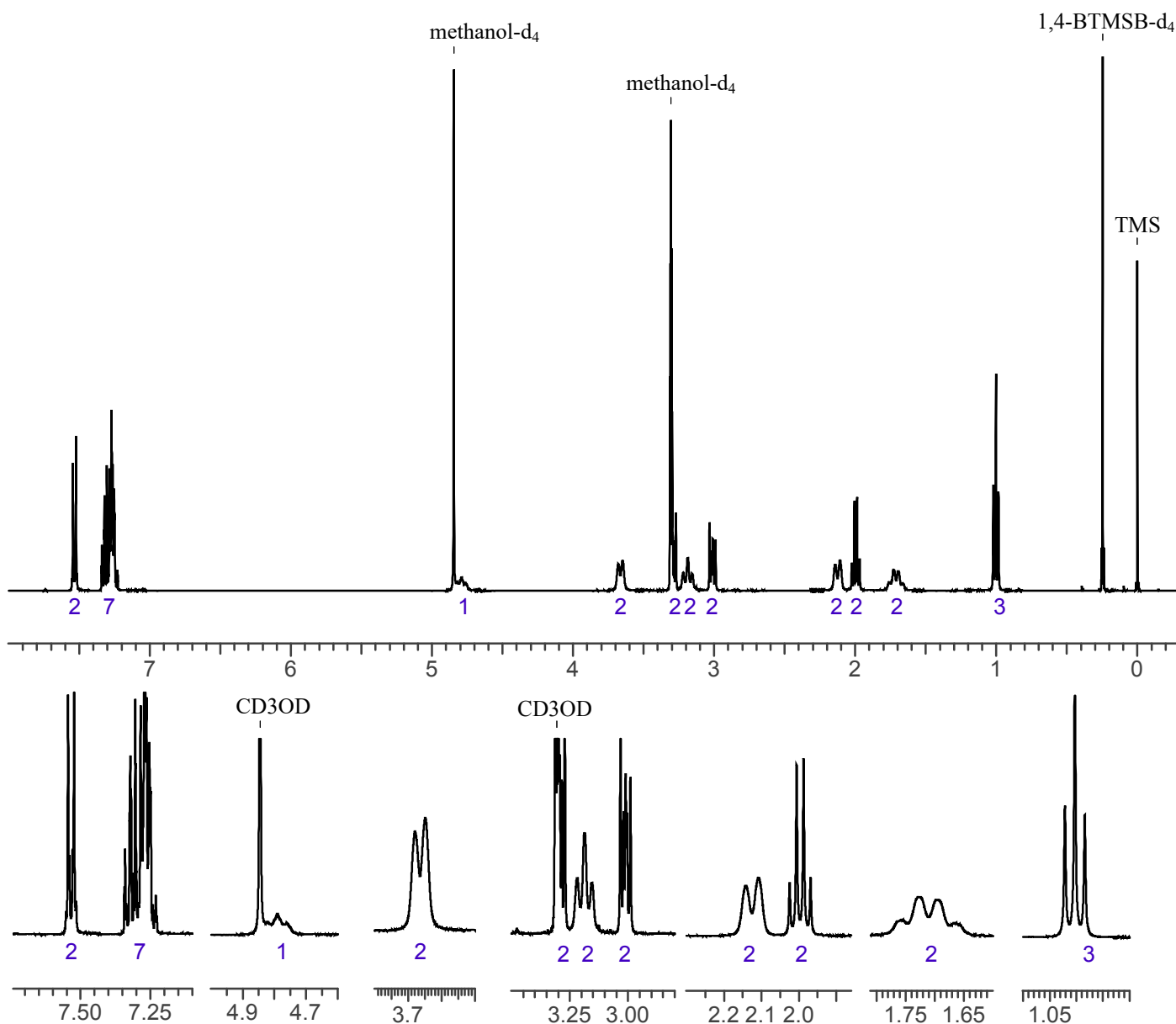



para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

Sample Preparation: Dilute analyte ~4 mg/mL into methanol.

Instrument: Agilent gas chromatograph operated in split mode with MS detector
Column: HP-5; 30m x 0.25 mm x 0.25 μ m
Carrier Gas: Helium at 1.5mL/min
Temperatures: Injector: 280°C
MSD transfer line: 280°C
MS Source: 230°C
MS Quad: 150°C
Oven program:
1) 100°C initial temperature for 1.0 min
2) Ramp to 280°C at 12 °C/min
3) Hold final temperature for 9.0 min
Injection Parameters: Split Ratio = 25:1, 1 μ L injected
MS Parameters: Mass scan range: 30-550 amu
Threshold: 100
Tune file: stune.u
Acquisition mode: scan
Retention Time: 17.817min

EI Mass Spectrum: para-Chlorofentanyl HCl; Lot# 0560132-2

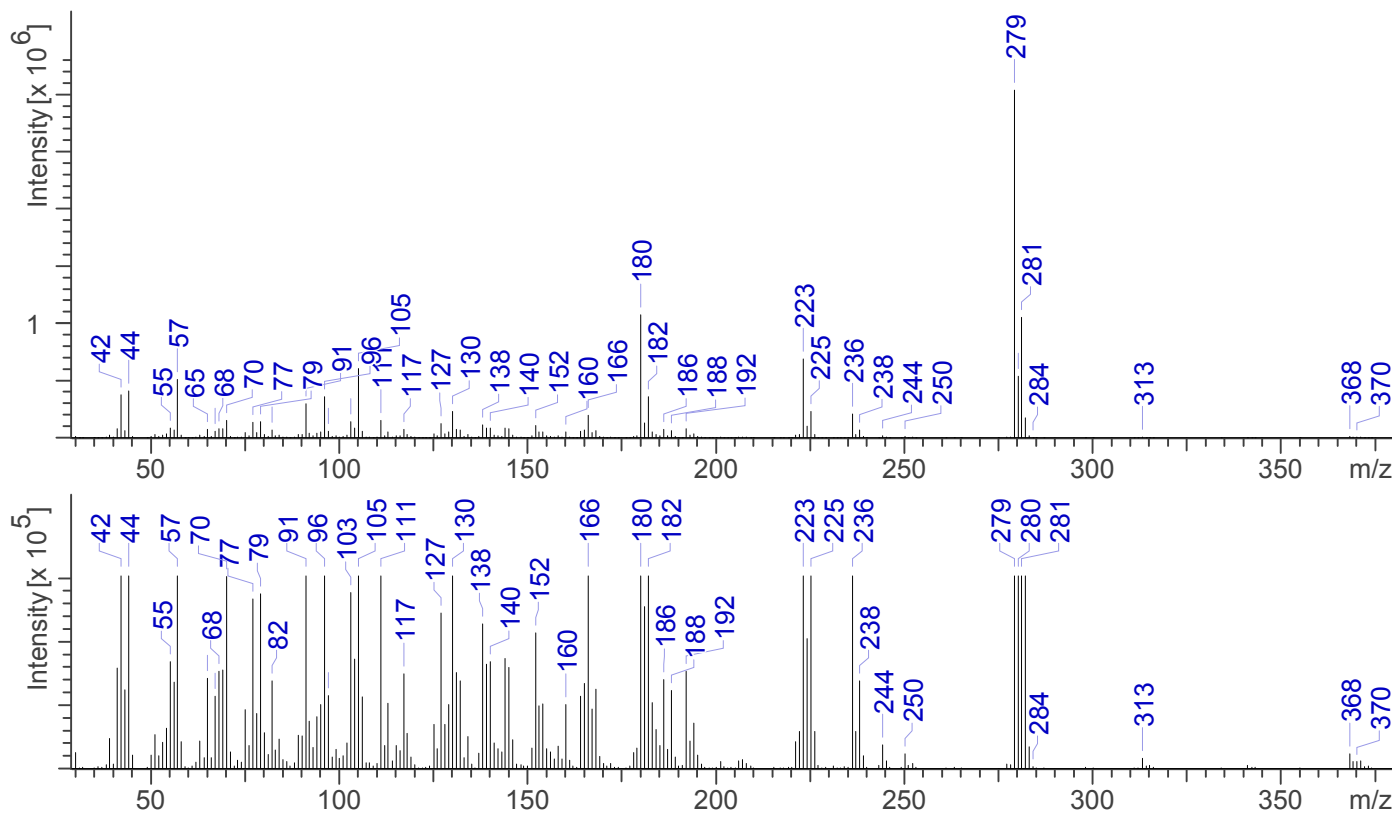



para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with diamond ATR attachment (3 bounce)

Scan Parameters:
Number of scans: 32
Number of background scans: 32
Resolution: 4 cm⁻¹
Sample gain: 4
Aperture: 150

FTIR-ATR (Diamond 3 Bounce): para-Chlorofentanyl HCl; Lot# 0560132-2

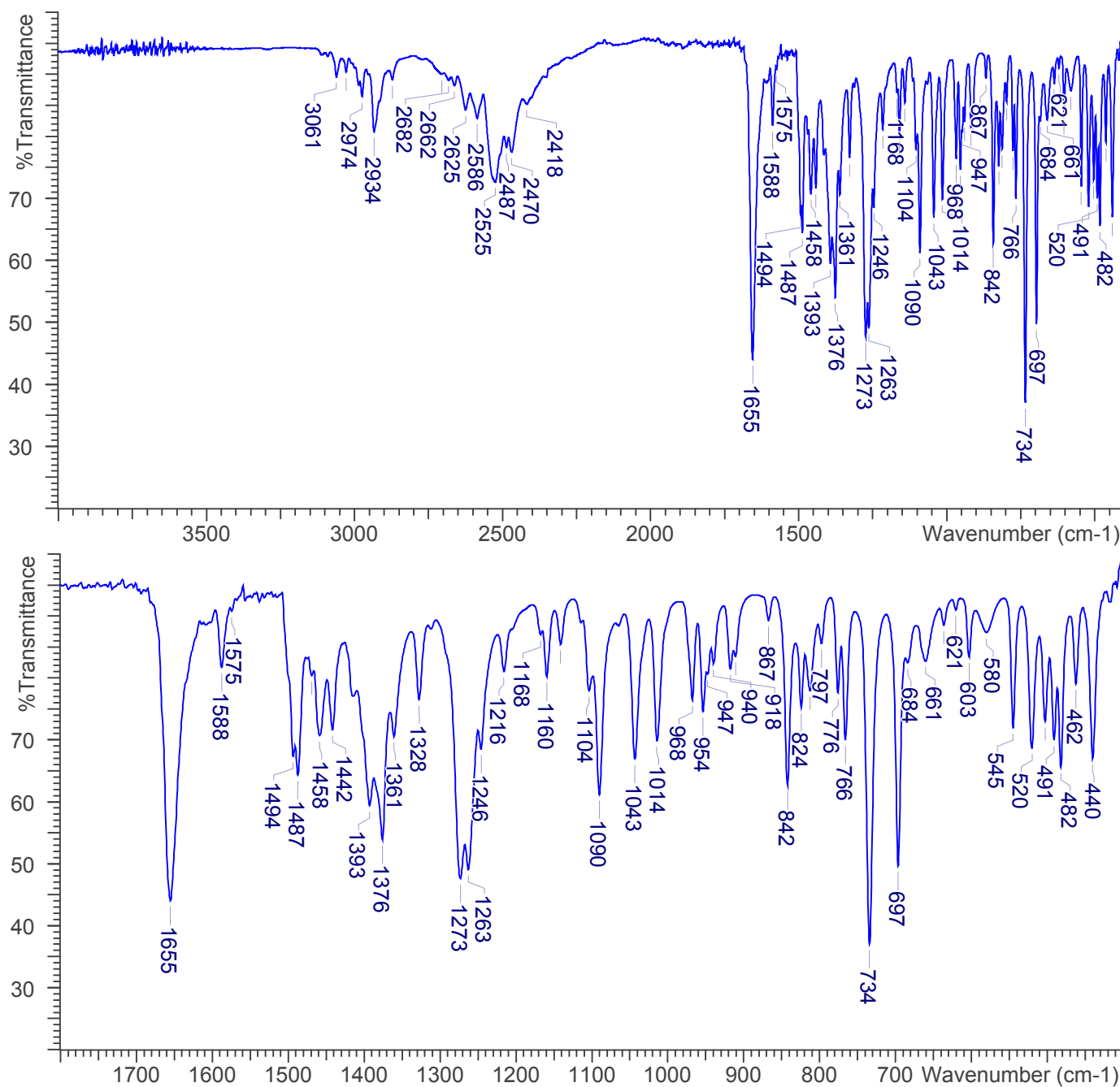



para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.





3.4. GAS CHROMATOGRAPHY /INFRARED DETECTION:

Sample Preparation: Dilute analyte 1mg/ml in CHCl_3

Instrument: Gas Chromatograph in split mode with Infrared Detection
Column: HP-5; 30m x 0.32 mm id x 0.25 μm
Carrier Gas: Helium at 2.0mL/min
Temperature: Injector: 280°C Split ratio: 2:1, 2 μl injection
65°C hold 0.5min., ramp to 310°C at 20°C/min., hold 5min.

IRD: Detector: Transfer line
Temp 280°C
Flow Cell Temp 280°C
Resolution 8 cm^{-1}

GC-IRD: para-Chlorofentanyl HCl; Lot 0560132-2

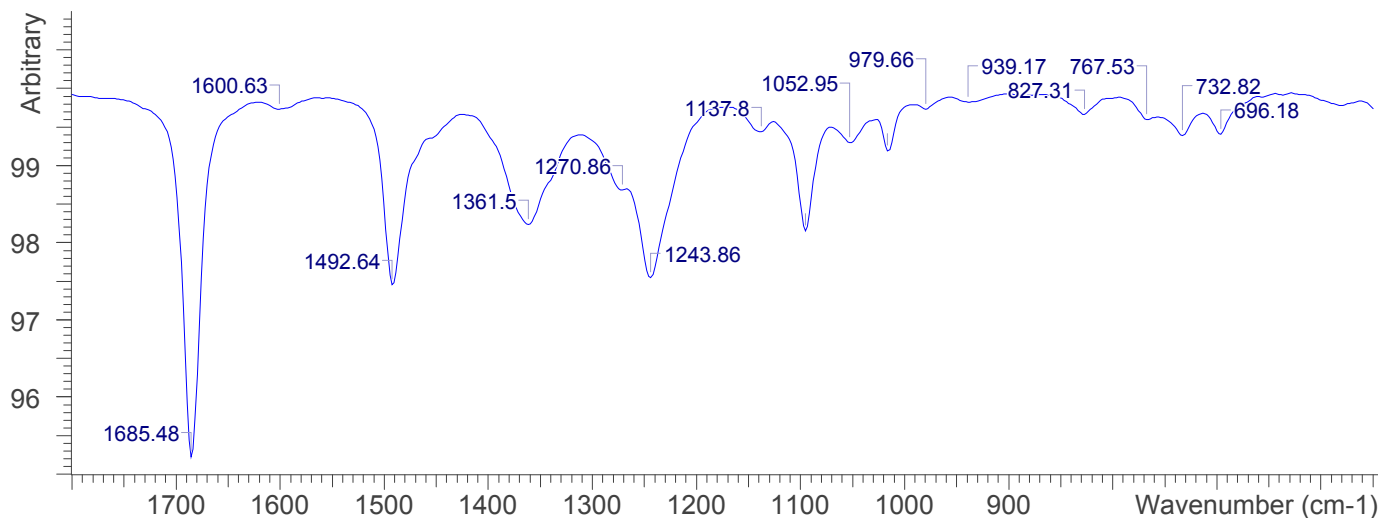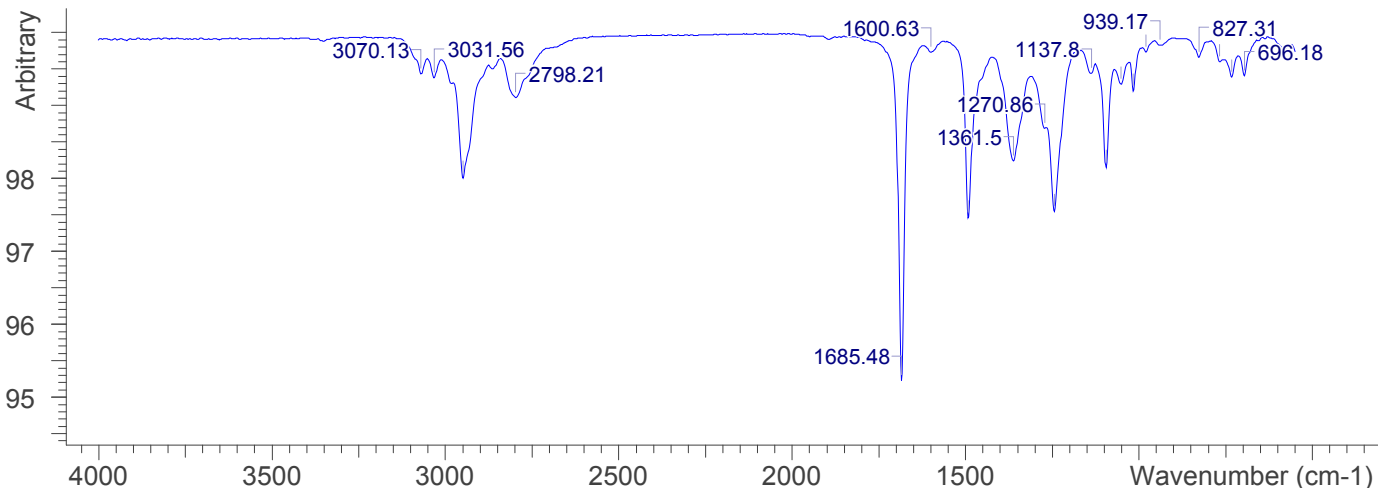



para-Chlorofentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



4. ADDITIONAL RESOURCES

No additional resources as of 07/01/2021